

Ultra-Care 4800



ALTERNATING PRESSURE/LOW-AIR LOSS MATTRESS SYSTEM

The **Ultra-Care 4800** is a combination mattress system providing both alternating pressure and low air loss for the management and treatment of pressure ulcers. Quick disconnect cells allow any area of the mattress to deflate and be zero pressure without affecting alternating benefits. In addition to the Alternating Pressure mode, the Static mode provides a firm surface when needed. Regardless of the operating mode, the mattress maintains the low air loss benefit.



GROUP 2

- Includes mattress, cover and pump
- 8" alternating pressure mattress with low air loss
- Pillow function maintains the air cells at the head for patient comfort
- Air cells made with nylon/polyurethane for increased durability and reliability
- Quick disconnect cells for easy cell replacement
- CPR quick deflate at head of mattress
- Quilted nylon zipper cover is low shear, fluid-resistant, and vapor permeable
- 8 liter-per-minute pump in an easy-to-use design
- Two modes of therapy - static or alternating
- Variable pressure setting for maximum comfort
- Support patients up to 350 lbs.
- HCPCS Code: E0277
- Two year warranty on mattress and pump

CELL CONSTRUCTION




- Twenty rows of 8" alternating pressure air cells, 9 with laser holes to keep patient cool and speed recovery
- Air cells made of nylon/polyurethane for increased durability and reliability
- Cells are easily cleaned

COVER



- Zipper cover is quilted vapor permeable nylon
- Water-resistant
- Easy cleaning
- Low shear/low friction

QUICK DISCONNECT CELLS

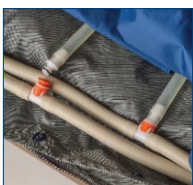

- Quick disconnects allow for easy cell replacement
- Quick disconnect cells allow any area of the mattress to be deflated and be zero pressure without affecting alternating benefits.

4581 PUMP



- 8 liter-per-minute pump
- Front panel switch provides two modes of therapy - static or alternating
- LED backlit switches with dial on panel for simple set up and operation
- Retractable hooks for bed use

Item #	Description	Wt. Cap.	Dimensions	Qty.
4800	Ultra-Care 4800 - 8" Mattress with 8 LPM Pump	350 lbs.	80" x 36" x 8"	each
4800-84	Ultra-Care 4800 - 8" Long Mattress with 8 LPM Pump	350 lbs.	84" x 36" x 8"	each
4800E	Ultra-Care 4800 - 8" Mattress with 16 LPM Pump (4253E)	450 lbs.	80" x 36" x 8"	each
4802 DC	Ultra-Care 4800 - 8" Mattress Only	350 lbs.	80" x 36" x 8"	each

Product Specifications:

ULTRA-CARE 4800 MATTRESS

Dimensions:	80" L 84" L x 36" H x 8" W
Cells:	20 rows of cells 8" high, 9 with laser holes for air loss benefit
Cell Connection:	Quick disconnect cells allow for single cell removal/replacement.
Inflation Time:	20 - 25 minutes with static inflate. 2 - 5 minutes with rapid inflate pump.
Material:	Nylon/Polyurethane for durability. Integrated foam base for additional protection.
Weight Capacity:	350 lbs. 450 lbs. with 4253E upgrade pump
Mattress Weight:	9 lbs.
Cover:	Urethane coated nylon that is vapor permeable and low shear/low friction
Base:	Reinforced nylon for durability. Bayonet type tie connectors to secure the mattress to the bed frame.
Transport Cap:	Located at end of the connector hose which when inserted, locks and prevents air loss from the mattress if the pump loses power or is disconnected.
Additional:	System uses modular components which allows for simple part replacement. CPR quick deflate at head of mattress.

ULTRA-CARE EXCEL 4581 PUMP

Input Rating:	110-120 VAC, 60Hz, 100mA, 8.5W
Power Cord:	SJT 14 Feet, 3 Prong
Air Output:	8 LPM
Indicator Lights:	Power On, Normal Pressure, Low Pressure, Mode Selection
Pressure Range:	20 - 70 mmHg \pm 10%
Cycle Time:	10 minutes
Dimensions:	11.8" L x 3.7" W x 7.5" H
Weight:	4.8 lbs.
Operating Modes:	Alternating Pressure or Static
Safety:	UL 544 & CSA

Optional Accessories:

16 LPM PUMP UPGRADE (ITEM #4253E)



- High-flow 16 liter-per-minute pump
- Dual compressor for increased patient comfort and weight capacity (100 lbs.)
- 10, 15 or 20 minute cycle times
- Five comfort settings
- Automatic reset from static to alternating after 45 minutes

SAFE-T-GUARD COVER (ITEM #SG8_FG)



- Perimeter cover reduces risk of falls and entrapment
- Foam borders are 4" high and sloped inward
- Zippered cover for secure fit
- Water resistant quilted top

QUICK INFLATE PUMP (ITEM #90501)



- Reduces inflation time dramatically
- 20 liter-per-minute flow pump
- Rapid deflate for quick transport
- Small and easy to use

* This is a Class II Medical Device. Federal law restricts these devices to sale by or on the order of a licensed health care professional.

* Providers are responsible for determining appropriate billing codes for Medicare Program claims. This information is included for your convenience and is not intended, nor should it be considered billing advice.

