

# Ultra-Care Xtra 4840

## BARIATRIC LOW AIR LOSS WITH CELL-ON-CELL MATTRESS DESIGN




The **Ultra-Care Xtra Bariatric Mattress System** is a combination mattress system providing both alternating pressure and low air loss for the for patients up to 800 lbs. This system features Stay-Flate™ cell-on-cell air bladder construction to provide a full ten inches of therapeutic support with bottom-out protection in the event of a power failure. In addition to the Alternating Pressure mode, the Static mode provides a firm surface when needed. Regardless of the operating mode, the mattress maintains the low air loss benefit.

**ULTRACARE**




- Includes mattress, cover and pump
- Stay-Flate™ cell-on-cell design provides air base and prevents "bottoming out" and up to 24 hours of power outage protection
- Pillow function maintains the air cells at the head for patient comfort
- Twenty rows of 10" alternating pressure air cells, nine with laser holes to keep the patient cool and speed recovery
- Air cells made with nylon/polyurethane for increased durability and reliability
- Quick disconnect cells for easy cell replacement
- CPR quick deflate at head of mattress
- Quilted nylon zipper cover is low shear, fluid-resistant, and vapor permeable
- High output digital pump
- Two modes of therapy - static or alternating
- Five pressure setting for maximum comfort
- Support patients up to 800 lbs.
- HCPCS Code: E0277
- Two year warranty on mattress and pump

### CELL CONSTRUCTION

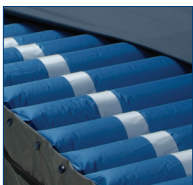

- Stay-Flate™ cell-on-cell design minimizes risk of bottoming out and provides air support during power outages
- Air cells made of nylon/polyurethane for increased durability and reliability
- Cells are easily cleaned

### QUICK DISCONNECT CELLS

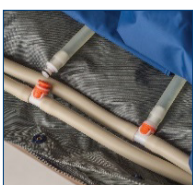

- Quick disconnects allow for easy cell replacement
- Quick disconnect cells allow any area of the mattress to be deflated and be zero pressure without affecting alternating benefits

### COVER



- Zipper cover is quilted vapor permeable nylon
- Water-resistant
- Easy cleaning
- Low shear/low friction

### 4253E PUMP



- High output digital pump
- Dual compressor for increased patient comfort and weight capacity
- 10, 15 or 20 minute cycle times
- Five comfort settings for easy pressure adjustment
- Automatic reset from static to alternating after 45 minutes

Item #	Description	Wt. Cap.	Dimensions	Qty.
4840-42	42" Ultra-Care Xtra System	800 lbs.	80" x 42" x 10"	each

## Product Specifications:

### ULTRA-CARE XTRA 4840 MATTRESS

Dimensions:	80" L x 42"W
Cells:	20 rows of cells 10" high, 9 with laser holes for air loss benefit. Seven inches of alternating air over three inches of static air. Cell-on-cell construction prevents "bottoming out" and up to 24 hours of power outage protection.
Cell Connection:	Quick disconnect cells allow for single cell removal/replacement.
Inflation Time:	25 - 30 minutes with static inflate. 2 - 5 minutes with rapid inflate pump.
Material:	Nylon/Polyurethane for durability
Weight Capacity:	800 lbs.
Mattress Weight:	25-35 lbs.
Cover:	Urethane coated nylon that is vapor permeable and low shear/low friction
Base:	Reinforced nylon for durability. Bayonet type tie connectors to secure the mattress to the bed frame.
Transport Cap:	Located at end of the connector hose which when inserted, locks and prevents air loss from the mattress if the pump loses power or is disconnected.
Additional:	System uses modular components which allows for simple part replacement. CPR quick deflate at head of mattress.

### 4253E DIGITAL PUMP

Input Rating:	AC 120V 60 Hz, 1 A
Power Cord:	SJT UL approved, 3-prong
Fuse Rating:	250V, 1A
Air Output:	12+ LPM
Pressure Range:	20 ~ 50 mmHg
Cycle Time:	10, 15, 20 cycle times selectable
Dimensions:	15.75" L x 5.5" W x 8" H
Weight:	7.9 lbs.
Operating Modes:	Alternating Pressure with low air loss or Static with low air loss. Automatic reset from static to alternating after 45 minutes.
Safety:	UL 544 & CSA
Additional:	Five comfort settings. Audible and visual alarms for air pressure and power supply. Retractable bed hooks for bed use.

## Optional Accessories:

### QUICK INFLATE PUMP (ITEM #90501)



- Reduces inflation time dramatically
- 20 liter-per-minute flow pump
- Rapid deflate for quick transport
- Small and easy to use

### SAFE-T-GUARD COVER (ITEM #SGU42\_FG)



- Perimeter cover reduces risk of falls and entrapment
- Foam borders are 4" high and sloped inward
- Zippered cover for secure fit
- Water resistant quilted top

\* This is a Class II Medical Device. Federal law restricts these devices to sale by or on the order of a licensed health care professional.

\* Providers are responsible for determining appropriate billing codes for Medicare Program claims. This information is included for your convenience and is not intended, nor should it be considered billing advice.