

Luna II CPAP & Auto CPAP Systems



The **3B Medical Luna II CPAP System** with RESlex technology and heated humidification is a single pressure system with advanced features and is a cost-effective option for treating sleep apnea. The Luna II Auto CPAP System uses unique Dual-Channel Quad-Phase technology. Increasing therapy comfort, RESlex Exhalation Relief eases the exhalation process making it easier to breathe out against the air flow. The Luna II humidifier chamber easily snaps into the side of the device using the same power supply as the machine. Optional Wi-Fi plug-ins offer the ability to capture patient compliance data from almost anywhere, anytime.

SPECIFICATIONS

| | |
|-----------------------|--|
| H2O Therapy Pressure: | 4-20 cm |
| Ramp Time: | 5-60 minutes <i>(5-minute increments)</i> |
| Dimensions: | 10.75" w x 7.25" d x 4.5" h <i>(device)</i> ; 14" x 14" x 6" <i>(shipping)</i> |
| Weight: | 4.2 lbs <i>(device)</i> ; 7.7 lbs <i>(shipping)</i> |
| Water Capacity: | 360 mL |
| Humidifier Settings: | Off, 1-5 (95°F-154°F / 35°C-68°C) |
| Humidifier Output: | ≥10 mg H2O/L |
| Warranty: | Two Years Limited <i>(device)</i> Three Months <i>(accessories)</i> |

FEATURES & BENEFITS

Large 3.5" LCD Display

Easily adjust settings and view last session to encourage user compliance.

Quiet Operation

with real time audible alert for tubing/mask leaks and automatic leakage/altitude compensation.

RESlex Exhalation Relief

with 3 comfort settings and inspiration trigger for automatic startup.

Remote Compliance Data

access with Wi-Fi or Cellular (optional), QR coding, iCode or SD card offering average statistics.

Embedded Memory


stores previous night and 365 night user raw data for up to 30 years!

iCode Connect

Register your patients to maintain records, generate sleep reports and manage HIPAA-compliance.



Travel Case

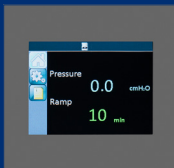

3.5" LCD Display

| Item # | Description | HCPCS | UOM | Case QTY |
|---------|--|--------------|------|----------|
| LG2000 | Luna II CPAP with Heated Humidifier | E0601, E0562 | Each | 1 |
| LG2A00 | Luna II Auto CPAP with Heated Humidifier | E0601, E0562 | Each | 1 |
| LG24510 | Luna II Replacement Humidifier Chamber | A7046 | Each | 1 |
| LG1020 | Luna II Reusable Foam Filter | A7038 | Pack | 6 |
| WL-100 | Luna WiFi Module | N/A | Each | 1 |
| WL-300 | Luna II 4G LTE Cellular Module | N/A | Each | 1 |