


PBSATINFOAM3680



ARUBA-200

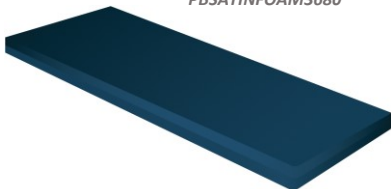

PBSMB Semi-Electric Bed Frame Assembly only

HBSMBEM ProBasics Semi-Electric Bed Ends w/Casters
 HBSMSS ProBasics Semi-Electric Single Motor, Head Spring, Foot Spring & Hi-Low Rod

PBSMB-FR Semi-Electric Bed w/ Full Rails (No Support Surface)

HBSMBEM ProBasics Semi-Electric Bed Ends w/Casters
 HBSMSS ProBasics Semi-Electric Single Motor, Head Spring, Foot Spring & Hi-Low Rod
 PB6034 ProBasics Full Bed Rails

PBSMB-HR Semi-Electric Bed w/ Half Rails (No Support Surface)

HBSMBEM ProBasics Semi-Electric Bed Ends w/Casters
 HBSMSS ProBasics Semi-Electric Single Motor, Head Spring, Foot Spring & Hi-Low Rod
 PB7035 ProBasics Half Bed Rails

PBSMB-FRFCPKG Semi-Electric Bed w/ Bed Ends, Full Rails & Fibercore Mattress

HBSMBEM ProBasics Semi-Electric Bed Ends w/Casters
 HBSMSS ProBasics Semi-Electric Single Motor, Head Spring, Foot Spring & Hi-Low Rod
 PB6034 ProBasics Full Bed Rails
 PBHBMATFB Fibercore Mattress

PBSMB-HRFCPKG Semi-Electric Bed w/ Bed Ends, Half Rails & Fibercore Mattress

HBSMBEM ProBasics Semi-Electric Bed Ends w/Casters
 HBSMSS ProBasics Semi-Electric Single Motor, Head Spring, Foot Spring & Hi-Low Rod
 PB7035 ProBasics Half Bed Rails
 PBHBMATFB Fibercore Mattress

PBSMB-FRISPKG Semi-Electric Bed w/ ProBasics Bed Ends, Full Rails & Innerspring Mattress

HBSMBEM ProBasics Semi-Electric Bed Ends w/Casters
 HBSMSS ProBasics Semi-Electric Single Motor, Head Spring, Foot Spring & Hi-Low Rod
 PB6034 ProBasics Full Bed Rails
 PBHBMATIS Innerspring Mattress

PBSMB-HRISPKG Semi-Electric Bed w/ Bed Ends, Half Rails & Innerspring Mattress

HBSMBEM ProBasics Semi-Electric Bed Ends w/Casters
 HBSMSS ProBasics Semi-Electric Single Motor, Head Spring, Foot Spring & Hi-Low Rod
 PB7035 ProBasics Half Bed Rails
 PBHBMATIS Innerspring Mattress

PBSMB-FRGP1PKG Semi-Electric Bed w/ Bed Ends, Full Rails & Group 1 Boxed Mattress

HBSMBEM ProBasics Semi-Electric Bed Ends w/Casters
 HBSMSS ProBasics Semi-Electric Single Motor, Head Spring, Foot Spring & Hi-Low Rod
 PB6034 ProBasics Full Bed Rails
 PBSATINFOAM3680 Group 1 Compressed Foam Mattress - Boxed

PBSMB-HRGP1PKG Semi-Electric Bed w/ Bed Ends, Half Rails & Group 1 Boxed Mattress

HBSMBEM ProBasics Semi-Electric Bed Ends w/Casters
 HBSMSS ProBasics Semi-Electric Single Motor, Head Spring, Foot Spring & Hi-Low Rod
 PB7035 ProBasics Half Bed Rails
 PBSATINFOAM3680 Group 1 Compressed Foam Mattress - Boxed

PBSMB-FRARBPKG Semi-Electric Bed w/ Bed Ends, Full Rails & Aruba 200 Mattress

HBSMBEM ProBasics Semi-Electric Bed Ends w/Casters
 HBSMSS ProBasics Semi-Electric Single Motor, Head Spring, Foot Spring & Hi-Low Rod
 PB6034 ProBasics Full Bed Rails
 ARUBA-200 ARUBA Group 1 Sculptured Foam Mattress - Boxed

PBSMB-HRARBPKG Semi-Electric Bed w/ Bed Ends, Half Rails & Aruba 200 Mattress

HBSMBEM ProBasics Semi-Electric Bed Ends w/Casters
 HBSMSS ProBasics Semi-Electric Single Motor, Head Spring, Foot Spring & Hi-Low Rod
 PB7035 ProBasics Half Bed Rails
 ARUBA-200 ARUBA Group 1 Sculptured Foam Mattress - Boxed



Lightweight Semi-Electric Bed

The ProBasics Lightweight Semi-Electric Bed is easy for one person to transport - with a foot section that weighs just 42 lbs. The robust, lightweight frame and snap on or off low-voltage DC motor box ensure quick, efficient set-ups and easy routine maintenance.

Your clients will take comfort that in the event of a power failure a 9V battery back-up may be used to allow movement to the lowest bed position. Plus, the smart, recessed design of the head and foot spring assemblies help reduce pinching and injuries.

FEATURES

- Easy to transport – foot section weighs just 42 lbs.
- Quick set-up & easy maintenance with motor & electronics that can be removed for service or cleaning
- Head and foot sections can be raised independently or simultaneously
- Large button, simple-to-use hand pendant
- Manual crank at the foot adjusts the frame height
- Scratch-resistant bed ends keep the bed looking new
- 9V battery back-up lets users go to the lowest bed position during power failures (batteries not included)
- Headspring support and foot assembly are recessed to help reduce pinching and injuries
- Trendelenberg & reverse Trendelenberg position capable
- Beds include 2 locking and 2 non-locking casters
- Quiet, efficient ETL-approved DC low voltage single motor box
- HCPCS Code: E0295***

PRODUCT SPECIFICATIONS

Surface Height:	15" min. – 20" max. (without mattress)	
Sleep Surface:	80" L x 36" W	
Overall Bed:	87" L x 36" W	
Weight Capacity:	450 lb*	
Limited Warranty:	5 year on frame** 2 year on motor - mechanical/electrical	

SHIPPING SPECIFICATIONS

Component:	Carton Size:	Shipping Weight:
Foot/Head Springs, Motor Box, pendant	45.5" x 10.5" x 37"	104 lbs.
Head/Foot Bed ends, Casters	38" x 7.5" x 32.25"	56 lbs.

ACCESSORIES

Works with ProBasics full/half rails, patient trapeze, bed extender and other components.

Item #	Description	UOM
PBSMB	ProBasics Semi-Electric Bed	each

* - Total distributed weigh capacity including user, mattress, bedding and any other items on the bed surface. Please refer to the Owner's Manual for complete details.

** - Limited lifetime warranty is offered on welds, please see Owner's Manual for complete warranty details.

*** - Providers are responsible for determining appropriate billing codes for Medicare Program claims. This information is only included for your convenience and is not intended, nor should be considered as Medicare billing advice.